


ADDITIONAL INFORMATION

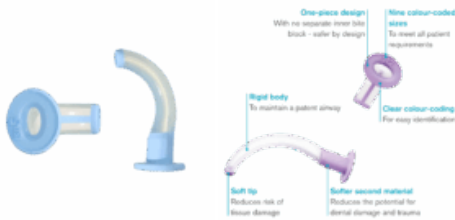

Weight	0.005 kg
Dimensions	11 × 6 × 2 cm
Carton Qty	25
Pricing Unit	Unit/1
Min Order Qty	1
ARTG ID	257365
GMDN Code	42424
Medical Device Class	Class 1s
UNSPSC Code	42271913
Tax Schedule	Non-Taxable
Is Medical Device	Yes
Is Serialized Item	No
Is Dangerous Good	No

PRODUCT DESCRIPTION

AeroRescue Guedel Airway Sterile is a one-piece oropharyngeal airway designed to help maintain an open airway in a patient. It is intended for use in emergency response, pre-hospital care, and clinical settings where rapid airway support equipment is required. The curved profile is formed to sit comfortably in the oral cavity while helping keep the tongue from obstructing airflow, supporting effective airway management when used by appropriately trained personnel.

Reliability in an airway adjunct starts with structural integrity. AeroRescue Guedel Airway Sterile uses a one-piece construction created through modern manufacturing methods that combine two materials with different physical characteristics into a single moulding. This design eliminates the separate inner bite block used in conventional Guedel airways. Removing this separate component helps reduce risks associated with loose or detached bite blocks, supporting safer handling and more consistent performance across demanding environments.

Comfort and patient safety are also addressed through material selection. A soft tip is incorporated to help reduce the risk of tissue damage during insertion. A softer second material is used in key areas to reduce the potential for dental damage and trauma, while the rigid body maintains the shape needed to help keep the airway patent. Together, these design elements support confident placement and stability when rapid airway access is needed.

Clear identification is essential when time is critical. AeroRescue Guedel Airway Sterile uses clear colour-coding to support quick selection and verification at a glance. The colour-coded sizing approach helps teams standardise equipment layout in resuscitation bags, airway roll kits, and treatment rooms, making it easier to locate the correct device quickly during assessments and patient

handovers.

Suitable for inclusion in rescue and airway kits, AeroRescue Guedel Airway Sterile provides a practical, purpose-designed solution for maintaining airway patency across a range of care environments, from first response to ongoing clinical support.

- **One-piece design:** eliminates a separate inner bite block for improved integrity.
- **Rigid body:** helps maintain a patent airway channel.
- **Soft tip:** designed to reduce the risk of tissue damage during insertion.
- **Softer second material:** helps reduce the potential for dental damage and trauma.
- **Clear colour-coding:** supports fast, simple size identification.

The **Size 00 (ISO 5.0)** 50mm option helps support appropriate patient fit, with clear colour-coding to make selection and verification faster when multiple sizes are kept in an airway kit.

This product is not a substitute for direct medical attention from a healthcare practitioner. Always read the label and follow the directions for use.

Blue Guedel Airway Size 00